

# Alaska Scientific Crime Detection Laboratory

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

### Document Structure

|   |    |
|---|----|
| <a href="#"><u>Introduction</u></a>   | 02 |
| <a href="#"><u>Abbreviations</u></a>  | 03 |
| <a href="#"><u>4.0 Management Requirements</u></a>                              | 06 |
| <a href="#"><u>5.0 Technical Requirements</u></a>                               | 19 |
| <a href="#"><u>Appendix A - - Physical Section Images Folder Flow Chart</u></a> | 26 |
| <a href="#"><u>Appendix B - - Processing Work Instructions</u></a>              | 27 |
| <a href="#"><u>Appendix C - - WIN/AAFIS, IAFIS</u></a>                          | 41 |
| <a href="#"><u>Appendix D - - DNA Sampling</u></a>                              | 42 |
| <a href="#"><u>Appendix E - - References</u></a>                                | 44 |
| <a href="#"><u>Appendix F - - Revision History</u></a>                          | 46 |

### Introduction

This document supplements the Alaska State Crime Laboratory Quality Assurance Manual. The guidelines and procedures in this manual are additional, Latent Print Discipline specific information.

The numbering scheme in this document follows that of the Alaska State Crime Laboratory Quality Assurance Manual. Supplemental requirements are found in Sections 4 and 5 of this document. Additional requirements in Sections 4 and 5 are listed by the Quality Assurance Manual criteria point they address. In Sections 4 and 5, the phrase "*Nothing additional*" means there is nothing additional to the requirements listed in the Laboratory Quality Assurance Manual.

ARCHIVED 7/27/2012

# Alaska Scientific Crime Detection Laboratory

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

### Abbreviations

|        |   |
|--------|---|
| 1°     | First Level Detail  |
| 2°     | Second Level Detail   |
| 3°     | Third Level Detail  |
| 10P    | Ten Prints or Ten Print Card  |
| =      | Control Negative  |
| +      | Control Positive  |
| #      | Number  |
| Ø      | Identification  |
| AAFIS  | Alaska Automated Fingerprint Identification System                  |
| AB     | Amido Black   |
| ACE-V  | Analysis, Comparison, Evaluation, Verification                      |
| AFIS   | Automated Fingerprint Identification System                         |
| ALS    | Alternate Light Source  |
| APIS   | Alaska Palm Identification System                                   |
| APSIN  | Alaska Public Safety Information Network                            |
| ACE-V  | Analysis, Comparison, Evaluation, Verification                      |
| B      | Blue  |
| BP     | Black Powder  |
| BTN    | Bearing the Name  |
| CA     | Cyanoacrylate Ester (Superglue)                                     |
| Cal    | Caliber   |
| CON    | Control   |
| CON=   | Control Negative  |
| CON+   | Control Positive  |
| CV     | Crystal Violet  |
| DAB    | Diaminobenzidine  |
| DET    | Detected  |
| DEV    | Developed   |
| DFO    | Diazafluoren-9-One  |
| DOB    | Date of Birth   |
| DL     | Driver's License  |
| ELIM   | Elimination   |
| ENV    | Envelope  |
| EVID   | Evidence  |
| ET     | Evidence Tape   |
| FLS    | Forensic Light Source   |
| FP     | Finger Print  |
| G      | Green   |
| GWL/OF | Green Wavelength Light with Orange Filter                           |
| HC     | Hard Copy   |
| IAFIS  | Integrated Automated Fingerprint Identification System (FBI system) |
| ID     | Identification  |
| IMP    | Impression  |
| IN     | Insufficient  |
| INC    | Inconclusive  |
| IND    | 1,2-Indanedione   |

# Alaska Scientific Crime Detection Laboratory

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012

Version: LP2012 R0

Effective: 01/24/2012

Archive: Archived

|       |   |
|-------|---|
| L#    | Lot Number  |
| LCV   | Leucocrystal Violet                               |
| LI    | Left Index Finger - Number 7                      |
| LL    | Left Little Finger - Number 10                    |
| LM    | Left Middle Finger - Number 8                     |
| LOV   | Latents of Value                                  |
| LP    | Latent Print                                      |
| LP/NV | Latents Present/No Value                          |
| LR    | Left Ring Finger - Number 9                       |
| LT    | Left Thumb Finger - Number 6                      |
| M:    | Marked  |
| MBD   | 7-pMethoxybenzylamino-4-nitrobenz-2oxa-1,3diazole |
| MENV  | Manila Envelope                                   |
| MP    | Magnetic Powder                                   |
| NA    | Nothing Additional                                |
| NAQ   | Not AFIS Quality                                  |
| NEG   | Negative  |
| NFE   | No Further Enhancement                            |
| NFRDD | No Further Ridge Detail Developed                 |
| NID   | No Identified / No Identification                 |
| NRD   | No Ridge Detail                                   |
| NSRDD | No Suitable Ridge Detail Developed                |
| NSRDP | No Suitable Ridge Detail Present                  |
| NV    | No Value  |
| NVRDD | No Visible Ridge Detail Developed                 |
| NVRDP | No Visible Ridge Detail Present                   |
| O     | Orange  |
| PD    | Physical Developer                                |
| POW   | Powder  |
| PREV  | Previous  |
| PROC  | Processed, processing                             |
| Pu    | Purple  |
| R     | Red   |
| R6G   | Rhodamine 6-G                                     |
| RAM   | Rhodamine, Ardrex and MBD                         |
| RD    | Ridge Detail                                      |
| RET   | Retained  |
| REV   | Reversed  |
| RI    | Right Index Finger - Number 2                     |
| RL    | Right Little Finger - Number 5                    |
| RM    | Right Middle Finger - Number 3                    |
| RR    | Right Ring Finger - Number 4                      |
| RT    | Right Thumb Finger - Number 1                     |
| RXN   | Reaction  |
| S/N   | Serial Number                                     |
| SRDD  | Suitable Ridge Detail Developed                   |
| SRDP  | Suitable Ridge Detail Present                     |
| SSN   | Social Security Number                            |
| STDS  | Standards   |

# Alaska Scientific Crime Detection Laboratory

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012

Version: LP2012 R0

Effective: 01/24/2012

Archive: Archived

|      |                                |
|------|--------------------------------|
| SUFF | Sufficient                     |
| T    | Teal                           |
| T/S  | Tape Sealed                    |
| UNS  | Unsuitable                     |
| UV   | Ultraviolet Light              |
| V    | Value                          |
| VIS  | Visual Examination             |
| VRDD | Visible Ridge Detail Developed |
| VRDP | Visible Ridge Detail Present   |
| W    | White                          |
| w/   | With                           |
| WIN  | Western Identification Network |
| WMP  | Magnetic Powder - White        |
| WP   | White Powder                   |
| WWb  | Wetwop - Black                 |
| WWw  | Wetwop - White                 |
| Y    | Yellow                         |
| ZC   | Zinc Chloride                  |

ARCHIVED 7/2/2012

### 4 Management requirements

#### 4.1 Organization

*Nothing additional*

#### 4.2 Management System

*Nothing additional*


#### 4.3 Document Control

*Nothing additional*

#### 4.4 Review of requests, tenders and contracts

- 4.4.1 The Latent Print Discipline Supervisor or designee reviews Request for Laboratory Services Forms and assigns casework to discipline analysts using the following steps.

The three adhoc queries used for case assignments are found in the list of available queries in the JTRAX Adhoc Query window shown in the screen shot below.

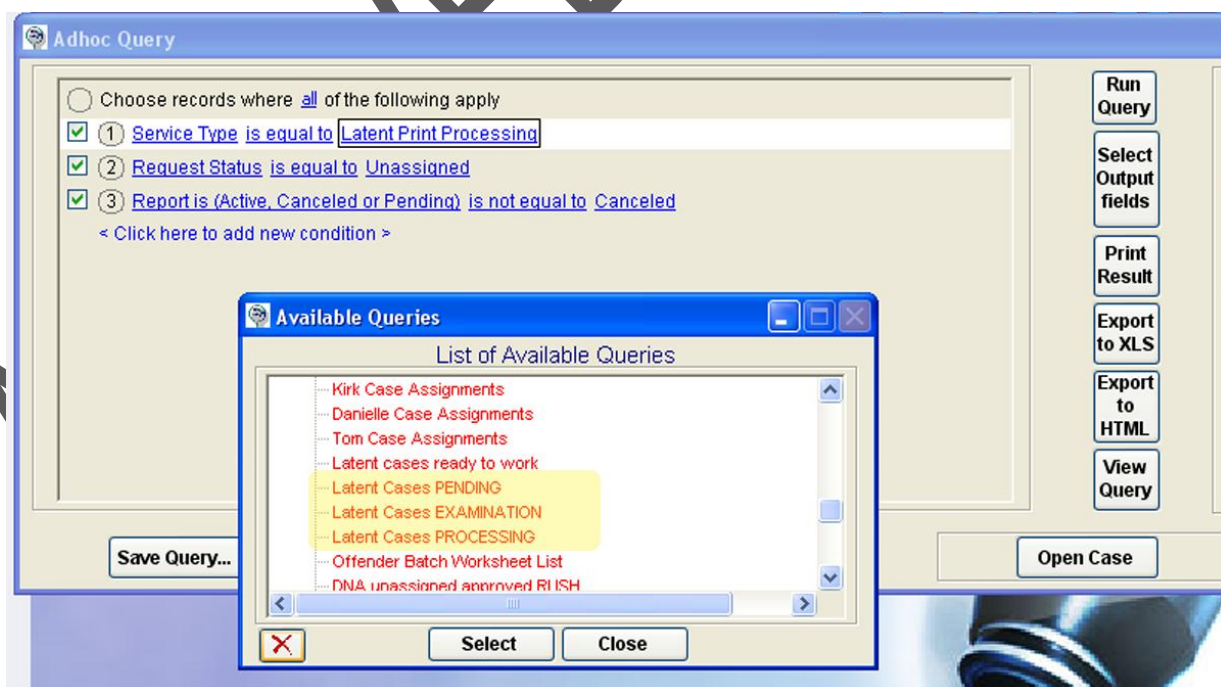

# Alaska Scientific Crime Detection Laboratory

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012

Version: LP2012 R0

Effective: 01/24/2012

Archive: Archived

The simplest request to assign is a Latent Print Examination request.

The **Latent Cases EXAMINATION** query will bring up all unassigned Latent Print Examination requests created after request items have been processed for latent prints. The majority of these unassigned Latent Print Examination requests are created by a Laboratory Technician who processed the request items. Latent Print Examination requests are ready for assignment to a Latent Analyst.

### LATENT CASES PENDING QUERY

The **Latent Cases PENDING** query will bring up all Pending Latent Print requests.

Latents Pending requests are created by evidence staff.

Open the RLS associated with each Latents Pending request.

If Latent Print analysis is the only type of analysis listed for case items, the request can be assigned using the following guidelines:

If the RLS indicates the items need to be processed for latent prints, the Latents Pending request is flipped to a Latent Print Processing request if assigned to a Laboratory Technician or the Latents Pending request is flipped to a Latent Examination request if assigned to a Latent Analyst.

If the RLS indicates there is no processing required (e.g. the items submitted are Lift card or photographic evidence), the Latents Pending request is flipped to a Latent Print Examination request and assigned to an Analyst in the Latent Print discipline. If the case only has lift cards, "**Lift Cards**" should be selected under **Reason** in the request window in JTRAX. This is shown in the screen shot below.

### REQUEST WINDOW

Request #0001 - Latent Print Examination

**Requesting Party Information**

Agency: ALASKA DEPARTMENT OF PUBLIC SAFETY Request Date: 09/02/2010

Badge Rep: Branch, Boyd

**Request Information**

Lab: LAB Analyst: Pippin, Turner

Department: Latents Due Date: / /

Service: Latent Print Examination Reason: Lift Cards

Requests: Complexity: APD Re-work, Approved Rush, Comm. Approved, LP/Bev ALC, Lift Cards, Rush

**Notes**

Requester: Assignor: Reviewer:

Related Evidence: 0 OK Cancel SOP

# Alaska Scientific Crime Detection Laboratory

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

Drugs, Arson and Biological Screening process items before Latent Prints. If Drugs, Arson or Biological Screening is requested for items that also have latent print work requested, the case is not assigned to a Latent Print Analyst until the other disciplines have completed their work.

### DRUG or ARSON/LATENT CASES

When the Drug discipline is finished working a case, the Drug Analyst will “flip” the Latents Pending request to a Latent Print Processing request. The **Latent Cases PROCESSING** query will bring up all unassigned Latent Print Processing requests created by the Drug Discipline. Before these requests are assigned, the Drug report should be reviewed and a note of which items were worked by the Drug discipline should be made in the **Assignor** notes box on the request window. This notes box is shown in the screenshot below.

Request #0001 - Latent Print Examination

**Requesting Party Information**

Agency: ALASKA DEPARTMENT OF PUBLIC SAFETY Request Date: 09/02/2010

Badge Rep: Branch, Boyd

**Request Information**

Lab: LAB Analyst: Pippin, Turner

Department: Latents Due Date: //

Service: Latent Print Examination Reason: Lift Cards

Requests: Complexity: APD Re-work, Approved Rush, Comm. Approved, LP/Bev ALC, Lift Cards, Rush

**Notes**

Requester: Assignor: Reviewer:

Related Evidence: 0 OK Cancel SOP

The Assignor notes box is also where any notes to the Laboratory Technician or Latent Print Analyst are entered. These notes could include a court date information, communication from the submitting agent, etc. Any information that the Case Assignor thinks the Assignee will need should be entered in the Assignor notes box.



## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

### BIOLOGICAL SCREENING/LATENT CASES

If a case has Biological Screening and Latent Print Processing assigned, review the RLS to see what type of screening is requested. If only Contact DNA swabbing and Latent Print processing is requested, the case is ready to be flipped to a Latent Examination request and assigned to a Latent Print Analyst. If Contact DNA swabbing AND other Biological Screening is requested for the case, check with the Supervisor of the Biology Discipline to clarify which discipline will do the contact DNA swabbing. A note should be put in the Assignor Notes box advising who will do the Contact DNA swabbing (e.g. "Latent Analyst to swab items 3, 4 and 6 prior to Latent Processing" or "Screening Discipline will swab for contact DNA").

Again, do not assign a case to a Latent Analyst until Biological Screening has completed their work on Biological Screening/Latent Processing evidence items.

4.4.2 *Nothing additional*

4.4.3 *Nothing additional*

4.4.4 *Nothing additional*

4.4.5 *Nothing additional*

### 4.5 Subcontracting of tests and calibrations

*Nothing additional*

### 4.6 Purchasing services and supplies

4.6.1 *Nothing additional*

4.6.2 Initial control testing of Rhodamine, Ninhydrin and other chemicals mixed at the laboratory is noted in the CHEM INV Excel Spreadsheet for each batch.

4.6.3 Purchasing documents will be saved online in the Uncontrolled Documents, Latent\_Share, Orders folder.

4.6.4 The Latent Discipline does not have any critical consumables.

### **4.7 Service to the Customer**

*Nothing additional*

### **4.8 Complaints**

*Nothing additional*

### **4.9 Control of nonconforming testing and/or calibration work**

*Nothing additional*

### **4.10 Improvement**

*Nothing additional*

### **4.11 Corrective Action**

*Nothing additional*

### **4.12 Preventive Action**

*Nothing additional*

ARCHIVED 7/2/2012

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

### 4.13 Control of Records

- 4.13.1.1 Latent print notes are printed to the JTRAX printer and then imported into JRTAX as a tif file.
- 4.13.1.2 Digital photos and/or scans for Latent Print casework are stored online in the PHYSICAL SECTION IMAGES folder. Inside this folder are two folders: Evidentiary Images and Work Product Images. Original images for cases are stored in the Evidentiary Images folder. Enhanced images and images to note Identifications, Comparisons and Verifications for cases are stored in the Work Product folder. Other documents such as scans of APSIN records, scans of ten print cards, etc. for cases are stored in the Work Product folder. A flow chart for the PHYSICAL SECTION IMAGES folder contents is stored in the PHYSICAL SECTION IMAGES folder. A copy of this flow chart is shown in **Appendix A**.
- 4.13.1.3 Access to the PHYSICAL SECTION IMAGES FOLDER and sub folders is limited. Latent Discipline Analysts, Crime Scene Technicians, the Latent Discipline Supervisor and the Crime Scene Technician Supervisor have access to the Evidentiary Images Folder. The Latent Discipline Supervisor and Latent Discipline Analysts have access to the Work Product Images folder.
- 4.13.1.4 *Nothing additional*
- 4.13.2.1 See 4.13.2.5.1 in this manual
- 4.13.2.2 See 4.13.2.5.1 in this manual
- 4.13.2.2.1 See 4.13.2.5.1 in this manual
- 4.13.2.3 *Nothing additional*
- 4.13.2.3.1 *Nothing additional*
- 4.13.2.3.2 *Nothing additional*
- 4.13.2.4 See 4.13.2.5.1 in this manual
- 4.13.2.5 *Nothing additional*

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

- 4.13.2.5.1 (Technical records in the Latent Print Discipline will meet the criteria as described in Appendix C Latent Print Examination Records) of ASCLD/LAB-*International* Supplement Standards.)

### **Start Dates – End Dates – Casework Activities**

The start date for casework is listed at the top of the first page of the case notes. The end date for casework is the date listed under “Conclusions” in the case notes. Dates for each process are noted next to the process. Processes occurring on the same date are noted in sequence from top to bottom or from left to right in an analyst’s notes.

### **Controls**

Any Positive or Negative control results for ninhydrin, DFO, IND, Cyanolacrylate or dyestain are documented next to the process in the analyst’s notes.

### **Analysis, Comparison, Evaluation, Verification**

The following methodology is a structured and systematic guide for comparing friction ridge detail. There are four parts to friction ridge identification methodology that includes: Analysis, Comparison, Evaluation, and Verification. This process, which is referred to as ACE-V, is repeated for each latent print developed. The primary purpose of these procedures is to establish unifying documentation for the methodology used in the comparison of friction ridge detail in the Latent Print Discipline.

The procedures presented are intended to assist the examiner in the comparison of friction ridge detail. They are to be used in conjunction with all applicable laboratory policies and proper scientific methodology.

For the following ACE-V procedure, a latent print is defined as friction ridge detail from an unidentified individual and a known exemplar is defined as friction ridge detail taken in a controlled manner from a known source, i.e. a “known” individual.

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

### Analysis

The suitability of latent prints is determined by analyzing three levels of friction ridge detail. This analysis considers the quality and quantity of the three levels of friction ridge detail.

Level 1 Detail (ridge flow) is not sufficient for individualization or exclusion. This level may include: general ridge flow, pattern configuration, core and delta location, distinction of finger versus palm, and other information enabling orientation.

Level 2 Detail (individual ridge path) enables individualization. This level may include ridge endings, bifurcations, dots, or combinations thereof.

Level 3 Detail (ridge shape) may enable individualization. This level may include: ridge width and shape, pores, edge contour, incipient ridges, breaks, creases, scars, etc.

If friction ridge characteristics are insufficient, the latent print is not suitable for identification purposes. Analysis is complete for that latent print.

If there are sufficient friction ridge characteristics, the latent print is determined to be suitable for identification purposes. The Analyst moves on to Comparison.

### Comparison

Conduct a comparison of the latent print to a known exemplar to determine if the ridge formations are in agreement. (Note: Comparisons can also be made between two latent prints or two known exemplars to determine if the prints came from the same source, i.e. individual. In these cases, one of the known exemplars or latent prints is chosen for the initial analysis and the process is identical to a latent print to known exemplar comparison).

If the friction ridge characteristics are not in agreement, the exemplar print is excluded as a source of the latent print, and the comparison process proceeds to other known exemplars for the case.

### Evaluation

If there is insufficient detail to form a conclusive identification.

1. The latent print may be compared to other prints, or
2. Determined to be inconclusive with documented justification, or
3. Determined to be not suitable for identification purposes with documented justification. In this instance analysis is discontinued.

If the prints have been determined to be in agreement and identified, the identification is documented on a composite of the latent print and the known exemplar. The friction ridge detail observed in agreement that support identification is marked on this composite. This composite is stored online in the Work Product Folder. An unmarked composite for verification is also stored in the Work Product Folder. At this point a "Latent Verification" request is created in JTRAX and assigned to another analyst.

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

### Exclusion of Individuals

Exclusion of an individual is the result of the comparison of two friction ridge impressions containing sufficient quality (clarity) and quantity of friction ridge detail which is not in agreement. Exclusion occurs when a latent print examiner, trained to competency, determines that two friction ridge impressions originated from different sources.

Two things are needed to exclude an individual.

1. An "anchor point" must be present to reliably determine the anatomical location of the latent print for an exclusion. An anchor point may include the following:
  - Delta
  - Core
  - Anatomical aspect allowing determination of origin location (i.e. outline of hand or finger, characteristic ridge flow or pattern)
  - Large field of ridge detail which may not have the above (i.e. hypothenar area of palm)
2. Clear known exemplar(s) from an individual that record ALL ridge detail that includes the "anchor point" present in the latent print.

All exclusions of individuals must be verified by a second analyst. All exemplars for an exclusion of an individual must be included in the Work Product folder with the latent print(s) excluded to an individual. A verification request is created for the exclusion of an individual.

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

### Verification

For identified latent prints, the assigned verifying analyst performs Analysis, Comparison, and Evaluation using the unmarked composite previously saved in the Work Product folder. The verifying analyst marks the friction ridge detail observed in agreement that support identification. This marked verification composite is stored online in the Work Product Folder.

For exclusions of individuals, the verifying analyst reviews the exemplars used by the original analyst to exclude the individual and performs Analysis, Comparison, and Evaluation to determine if they agree with the exclusion of the individual.

If the verifying analyst does not agree with an identification or exclusion of an individual, then the Latent Print Supervisor is notified for resolution. The Latent Print Supervisor may designate another analyst to review the latent in question.

Additional reviews are performed until such a time as the identification or exclusion is confirmed or refuted.

If identification is refuted, appropriate corrective actions are initiated by the Latent Discipline Supervisor.

### **Image Security**

Access to the PHYSICAL SECTION IMAGES FOLDERS is limited as noted previously in this manual under 4.13.1.3

4.13.2.5.2 There are no instrumental analyses in the Latent Print Discipline.

4.13.2.6 *Nothing additional*

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

- 4.13.2.7 Typical casework reporting should follow the recommended reporting statements as appropriate. Not every situation can be represented by the listed report examples, so report statements may include but are not limited to the examples given.

### Latent Processing Report Guidelines

1. List description of all items on which processing was performed.
2. Latent Discipline Laboratory Technicians indicate if areas of ridge detail were developed and documented.
3. Latent Discipline Analysts indicate if latent prints were developed and suitable for identification purposes.
4. Indicate disposition of retained items, submitted exemplars, and any preserved latent prints (e.g. digital images, lifts, etc.).

Examples of Latent Processing Statements include but are not limited to the following:

"No latent prints were developed on [item]." (No ridge detail)

"No suitable latent prints were developed on [item]." (Ridge detail present not suitable for identification)

"Suitable latent prints were developed and preserved (photographed, scanned, or lifted) on [item]." (Ridge detail present suitable for identification)

"The suitable latent print developed appears to be a palm print (or finger, tip, side, joint, palm, toe, sole, etc)."

"Please submit known (fingerprints, palm prints or foot prints) of individuals associated with the case for comparisons to be made."



## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

### Latent Print Comparison (Identification/Exclusion/Inconclusive) Report Guidelines

Reports regarding latent print comparisons must include:

1. Item description **[item]** (include item number, item description, and latent location as necessary).
2. Finger/palm/foot identified **[print designation]**.
3. Source of Exemplar **[subject name]** (include subject name and associated APSIN ID#, DOB, or other identifier as necessary).

#### Examples of Identification Statements include but are not limited to the following:

"The latent print present on **[item]** matches the **[print designation]** (fingerprint, palm print, toe print, sole print) on the card (submitted or on file at SCDL) bearing the name of **[subject name]**."

"The latent print present on **[item]** has been identified to the **[print designation]** (fingerprint, palm print, toe print, sole print) on the known prints (submitted or on file at SCDL) bearing the name of **[subject name]**."

#### Examples for Identifications Made as a Result of an Automated Search

"As a result of an AFIS (Automated Fingerprint Identification System) search, the latent has been identified as being the same as the **[print designation]** (fingerprint, palm print) on the known prints bearing the name of **[subject name]**."

"The latent print present on **[item]** has been identified to the **[print designation]** (fingerprint, palm print) on the known prints bearing the name of **[subject name]**, obtained during an automated search."

#### Examples of Exclusion to Exemplar Statements

"The latent print does not match the fingerprints (palm prints or foot prints) on the card bearing the name of **[subject name]**."

"The latent print present on **[item]** was not identified to the known prints (submitted or on file at DPS) bearing the name of **[subject name]**."

"The latent print was compared to known prints (submitted or on file at DPS) bearing the name of **[subject name]** with negative results."

#### Examples of Exclusion of Individual Statements

"**[subject name]** is excluded as a source for the latent print found on **[item]**."

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

### Examples of Inconclusive Statements

"The latent print cannot be identified or excluded as being made by any of the known prints on the card bearing the name of **[subject name]**, due to **[reason]**".

"The latent print present on **[item]** was not identified to the known prints (submitted or on file) bearing the name of **[subject name]**. However, the known prints do not fully represent the area of skin needed for a complete examination to be made. Submit a fully rolled set of (fingerprints/palm prints/to include joints of fingers, etc.) for a complete comparison to be made."

"The known prints DPS has on file or the known prints submitted are not suitable for a complete examination to be made."

### **AFIS Database Searches Report Guidelines**

#### Examples of AFIS Statements:

"An AFIS search on the remaining unidentified latent prints has been performed with negative results. The latent prints have been registered in the WIN/AAFIS database. In the event that a hit is generated at a later date, an additional report will follow."

"The suitable latent prints that (were not suitable for entry into AFIS/did not meet the criteria for an AFIS search) were not identified to known exemplars bearing the name of **[subject name]**."

4.13.2.8 *Nothing additional*

4.13.2.9 *Nothing additional*

4.13.2.10 *Nothing additional*

4.13.2.11 *Nothing additional*

4.13.2.12 All identifications and exclusions must be verified.

4.13.2.13 Latent Print Discipline abbreviations are listed in the Abbreviations Section of this manual.

### **4.14 Internal Audits**

*Nothing additional*

### **4.15 Management Reviews**

*Nothing additional*

### 5 Technical requirements

#### 5.1 General

- 5.1.1 Laboratory Technicians assigned work at the laboratory in the Latent Discipline and new Analysts are competent to process latent print evidence after successfully completing Section 1 of the Latent Print Training Manual. New Analysts are competent to process latent print evidence and compare latent prints after successfully completing Section 1 and Section 2 of the Latent Print Training Manual.
- 5.1.2 There are no measurements of uncertainty in the Latent Print Discipline
- 5.1.3 The Cyanolacrylate/Dyestain process is control tested each time it is used. A fingerprint is placed on a clear piece of glass or plastic and processed with evidence items. Ninhydrin, DFO and 1,2 IND processes are control tested each time they are used. A fingerprint is placed on a clean, white sheet of paper and processed with evidence items. Any Positive or Negative control results are recorded in the Analyst's case notes.
- 5.1.3.1 Initial control testing of Rhodamine, Ninhydrin and other chemicals mixed at the laboratory is noted in the CHEM INV Excel Spreadsheet for each batch.

#### 5.2 Personnel

- 5.2.1 All training documented in the Latent Print Training Manual is supervised by competent, experienced Forensic Scientist III or Forensic Scientist IV analysts
- 5.2.1.1 All training documented in the Latent Print Training Manual is signed off by competent, experienced Forensic Scientist III or Forensic Scientist IV analysts
- 5.2.1.2 The Latent Print Training Manual includes a section on court testimony. A moot court is required before an analyst is released for independent casework.
- 5.2.1.2 The Latent Print Training Manual includes sections on other forensic disciplines, court procedures and ethics.
- 5.2.2 *Nothing additional*
- 5.2.3 *Nothing additional*
- 5.2.4 *Nothing additional*
- 5.2.5 *Nothing additional*
- 5.2.6 *Nothing additional*
- 5.2.7 *Nothing additional*

#### 5.3 Accommodation and environmental conditions

*Nothing additional*

### 5.4 Test and calibration methods and method validation

5.4.1 The accepted processing methods used in the Latent Discipline are:

Amido Black  
Cyanoacrylate  
DFO  
IND  
Ninhydrin  
Physical Developer  
Powders - Plain (All colors and fluorescent powders from approved providers)  
Powders - Magnetic (All colors)  
Rhodamine  
Wetwop (white and black)

Mixing instructions for reagents made at the laboratory are in the CHEM INV Excel Spreadsheet.

Purchased reagents such as Cyanoacrylate, Physical Developer, Powders, Wetwop, are purchased from an approved vendor.

5.4.2 Processing used for a case evidence is left to analyst discretion. Processing guidelines are listed in **Appendix B** of this manual.

5.4.3 *Nothing Additional*

5.4.4 *Nothing Additional*

#### 5.4.5 Validation of Methods

5.4.5.1 *Nothing Additional*

5.4.5.2 *Validation records are stored in the Validations folder in the Latents Share folder.*

5.4.5.3 *Nothing Additional*

5.4.5.4 Performance Check records are stored in the Performance Checks folder in the Latents Share folder.

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

### 5.4.6 Estimation of uncertainty of measurement

5.4.6.1 *Nothing Additional*

5.4.6.2 Measurement of uncertainty does not apply to the Latent Print Discipline.

5.4.6.3 *Nothing Additional*

### 5.4.7 Control of data

5.4.7.1 *Nothing Additional*

5.4.7.2 *Nothing Additional*

5.4.7.2.1 Access to the PHYSICAL SECTION IMAGES FOLDER and sub folders is limited. Latent Discipline Analysts, Crime Scene Technicians, the Latent Discipline Supervisor and the Crime Scene Technician Supervisor have access to the Evidentiary Images Folder. The Latent Discipline Supervisor and Latent Discipline Analysts have access to the Work Product Images folder.

## 5.5 Equipment

5.5.1 Equipment used in the Latent Discipline consists of:

Ninhydrin Humidity Chamber

Cyanoacrylate Chamber

DFO Oven

532 nm Light Source

Digital Cameras

Balance

5.5.2 The Ninhydrin Humidity Chamber, Cyanoacrylate Chamber, 532 nm Light Source and DFO oven are control tested when evidence items are processed to ensure proper function.

5.5.3 Equipment manuals are stored in the Latents Share folder.

5.5.4 *Nothing Additional*

5.5.5 Spread sheets with equipment records are stored in the Latents Share folder.

5.5.6 Balances used for chemical preparation in the Latent Print Discipline are checked/calibrated yearly by an approved outside vendor. Normal maintenance includes keeping the balance clean and leveled.

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

5.5.7 *Nothing Additional*

5.5.8 The vendor that checks the chemical preparation balance attaches a sticker verifying the balance was checked.

5.5.9 *Nothing Additional*

5.5.10 *Nothing Additional*

5.5.11 *Nothing Additional*

5.5.12 *Nothing Additional*

### 5.6 Measurement traceability

#### 5.6.1 General

5.6.1.1 *Nothing Additional*

5.6.2 *Nothing Additional*

#### 5.6.3 Reference Standards and Reference Materials

5.6.3.1 *Nothing Additional*

5.6.3.2 *Nothing Additional*

5.6.3.2.1 The laboratory has a known palm reference collection. It consists of palm cards submitted for casework. Palm cards are identified by laboratory case number and name and filed alphabetically by last name. The palm collection is stored in the archived latent case file room which can only be accessed by Latent Discipline personnel. The key for the key for the latent case archive room is kept in a locked key box in the section.

5.6.3.3 *Nothing Additional*

5.6.3.4 *Nothing Additional*

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

### 5.7 Sampling

- 5.7.1 The question being asked in the Latent Print Discipline is "Did an individual leave a latent print on an item of evidence?" An analyst can select samples by quality and/or quantity. An analyst can select the evidence that has the best chance of retaining latent prints or an analyst can select and process a portion of the evidence until the individuals of interest are identified. A technician can employ sample selection if they consult with an analyst. Once an analyst identifies the individual(s) of interest on an item, processing and comparison can cease.
- 5.7.2 Unusual sample selection situations must be approved by the latent discipline supervisor.
- 5.7.3 An analyst or laboratory technician must document sampled selection in their notes. When sample selection is employed, the selected items must be identifiable at a later date. An analyst or laboratory technician must initial the sample selection items in a case at a minimum. If this is not possible, the sampled selection items may be repackaged and the analyst or technician must, at a minimum, initial the packaging.

### 5.8 Handling of test and calibration items

- 5.8.1 *Nothing Additional*
- 5.8.1.1 *Nothing Additional*
- 5.8.1.1.1 *Nothing Additional*
- 5.8.1.1.2 *Nothing Additional*
- 5.8.2 *Nothing Additional*
- 5.8.3 *Nothing Additional*
- 5.8.4 *Nothing Additional*
- 5.8.4.1 *Nothing Additional*
- 5.8.4.2 Evidence not in the process of examination is stored in a locked tote or evidence cabinet. Larger evidence may be stored in the latent case archive room.
- 5.8.4.3 *Nothing Additional*
- 5.8.4.4 *Nothing Additional*
- 5.8.4.5 *Nothing Additional*
- 5.8.4.6 Analysts using the WIN/AAFIS system will use AFIS Latent Fingerprint Best Practices, Western Identification Network, Inc., October 2008 as a reference. A copy of this guide is located in the WIN-AAFIS Folder in the Latent\_Share Folder.



## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

5.8.4.6.1 Palm files are treated as Reference Material

5.8.4.6.1.a *Nothing Additional*

5.8.4.6.1.b *Nothing Additional*

5.8.4.6.2 Palm cards are identified by laboratory case number and name

5.8.4.6.3 The palm collection is stored in the locked archived latent case file room. The key for the latent case archive room is kept in a locked key box in the section

5.8.4.6.3 The locked key box can only be accessed by Latent Discipline personnel.

### 5.9 Assuring the quality of test and calibration results

5.9.1 *Nothing Additional*

5.9.1.1 Control testing is covered under criteria point 5.1.3 in this manual.

5.9.2 If a control test is negative, the Latent Discipline supervisor will be notified. The Latent Discipline supervisor or a designee will take any necessary corrective action.

5.9.3 Each Latent Print analyst will take the CTS proficiency test yearly.

5.9.3.1 *Nothing Additional*

5.9.3.2 *Nothing Additional*

5.9.3.3 Each Latent Print analyst will take the CTS proficiency test yearly.

5.9.3.3.1 *Nothing Additional*

5.9.3.3.2 There is no external latent print processing proficiency test. Each Laboratory Technician processing evidence items and issuing reports will take an internal, latent print processing proficiency test approved and/or prepared by the Latent Discipline Supervisor.

5.9.3.4 Each Latent Print analyst will take the CTS proficiency test yearly.

5.9.3.5 *Nothing Additional*

5.9.3.6 *Nothing Additional*

5.9.4 *Nothing Additional*

5.9.4.1 In addition to the guidelines listed in the Laboratory Quality Assurance Manual, a Latent Print Discipline technical review will include:

Appropriateness of processes used based on substrate and possible composition of latents (e.g. bloody impressions visible indicating a blood reagent as a processing choice).

Proper process sequence

5.9.4.2 *Nothing Additional*

5.9.4.3 *Nothing Additional*



## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012  
Effective: 01/24/2012

Version: LP2012 R0  
Archive: Archived

5.9.5 *Nothing Additional*

5.9.5.1 In addition to the guidelines listed in the Laboratory Quality Assurance Manual, a Latent Print Discipline administrative review will include:

Check that the photos listed in the case notes match the photos stored in the PHYSICAL SECTION IMAGES folder for the case.

5.9.6 *Nothing Additional*

5.9.7 *Nothing Additional*

### 5.10 Reporting the results

5.10.1 *Nothing Additional*

5.10.2 *Nothing Additional*

#### 5.10.3 Test Results

5.10.3.1 *Nothing Additional*

5.10.3.2 *Nothing Additional*

5.10.3.3 *Nothing Additional*

5.10.3.4 *Nothing Additional*

5.10.3.5 Results and report wordings are covered under criteria point 4.13.2.7 in this manual.

5.10.3.6 Exclusions and report wordings are covered under criteria point 4.13.2.7 in this manual.

5.10.3.7 Inconclusive results and report wordings are covered under criteria point 4.13.2.7 in this manual.

5.10.4 *Nothing Additional*

5.10.5 The basis upon which opinions and interpretations are made is documented in the Identification and Verification composites stored in the Work Product folder in the PHYSICAL SECTION IMAGES folder and in the Latent Discipline analyst's notes for the case.

5.10.6 *Nothing Additional*

5.10.7 *Nothing Additional*

5.10.8 *Nothing Additional*

5.10.9 *Nothing Additional*

# Alaska Scientific Crime Detection Laboratory

## Latent Print Discipline - Additional Guidelines and Procedures

Issued: 01/23/2012


Effective: 01/24/2012

Version: LP2012 R0

Archive: Archived

### APPENDIX A - - PHYSICAL SECTION IMAGES FOLDER FLOW CHART

#### PHYSICAL SECTION IMAGE FOLDER OUTLINE

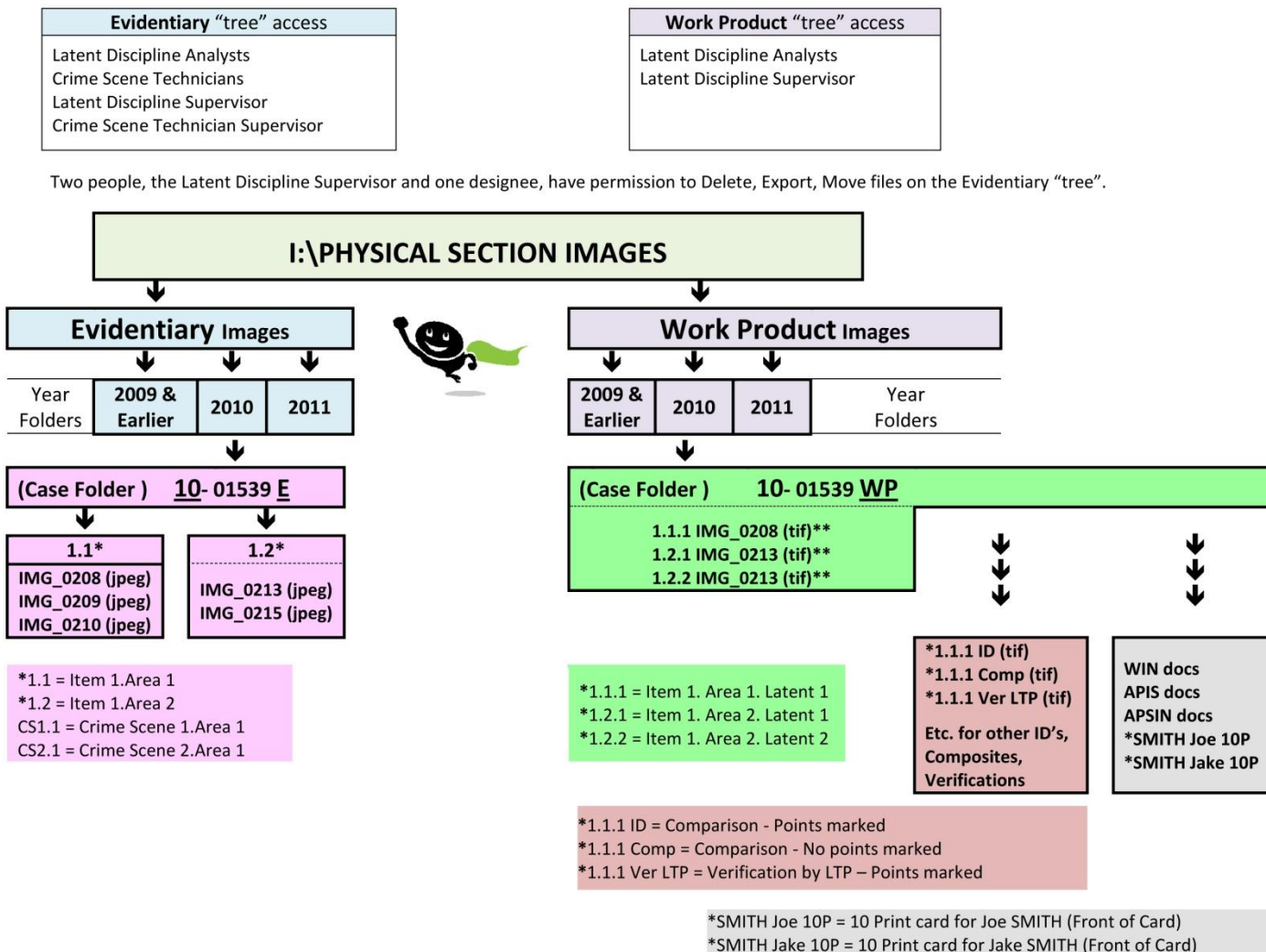

## APPENDIX B - - PROCESSING WORK INSTRUCTIONS

This appendix describes the use and safety for the different processing techniques used in the Latent Print Discipline.

Instructions for preparing chemicals made at the laboratory are found in the JTRAX Chemical Inventory. Exact measurements and proportions when preparing chemical solutions are desirable for consistent quality, but successful results in developing latent fingerprints are not dependent upon unequivocal accuracy. There is considerable latitude in preparing chemical solutions for latent fingerprint techniques without adversely affecting the successful development of latent prints.

ARCHIVED 7/2/2012

## **AMIDO BLACK**

### **Description of Process**

Amido Black, also known as naphthol blue-black, is used to develop or enhance latent prints that have been left in blood. Amido Black stains the proteins in the blood turning the print a dark blue or black color. This is not the only blood print development technique available.

Amido Black may destroy blood for serology/DNA testing. Have evidentiary blood samples chemically tested and preserved by appropriate personnel prior to processing. It will not develop prints in perspiration, fats and oils, or salts. The background of a porous item may also stain, causing weak bloody prints to not be detected.

Cyanoacrylate ester fuming may be detrimental to this process.

### **Safety**

Mix only in a vent hood.

When mixing or using, must wear gloves and eye protection.

Chemicals are flammable and skin irritant.

Caution should always be exercised around a bloody crime scene or handling items which contain blood.

Excess is disposed of as any flammable liquid.

Protective lab coats, footwear, eyewear, and latex gloves should be worn.

Since Amido Black is mixed with methanol, which is highly flammable, extreme caution should be taken when at a crime scene as to make sure.

## Process for Use

Amido Black is typically prepared in the laboratory and not purchased as a working solution. Shelf life is indefinite for Amido Black and Rinse Solution.

1. Preserve any suitable visible prints present on evidence prior to applying Amido Black solutions.
2. Apply the Amido Black base solution by dipping, spraying, or using a squirt bottle to dried prints in apparent blood. Apply until the entire print has turned from a reddish-brown color to a blue-black color. Background staining may occur.
3. If necessary, the base solution can be re-applied before the final rinse to achieve sufficient clarity.
4. Rinse off excess base solution with the rinse solution (use additional rinses as necessary to achieve sufficient clarity).
5. Let dry.
6. Examine the item for latent prints and indicate the latent with suitable markings as appropriate to be preserved by digital imaging/photography.

Note: Developed latent prints on some dark-colored surfaces may be viewed with a light source for increased contrast.

## Sequence

Amido Black is typically utilized instead of other processes (cyanoacrylate ester fuming, powders, etc.). Although a light application of cyanoacrylate ester fuming may be applied previous to Amido Black application to preserve latent prints not in apparent blood.

# CYANOACRYLATE

## Description of Process

The super glue process was first used for latent print development by the Japanese police in the late 1970s. A glue containing cyanoacrylate is placed in an airtight chamber with the evidence. As the glue dries, the fumes from the drying glue circulate throughout the chamber adhering to the latent print residue left on the evidence. The prints may then be dusted and lifted or preserved through photography and/or digital imaging. The process will develop fresh as well as old prints.

## Safety

Precautions should be taken as to not get the glue on your skin. Wear eye protection and latex gloves. If you do get glue on your skin and get attached to something, do not try to pull apart. Use water or acetone and then rub apart to release. Use in a vent hood or use an exhaust system to remove fumes from the chamber prior to opening the fuming chamber and removing the evidence. Cyanoacrylate ester fuming may be a respiratory and eye irritant.

## Process for Use

Cyanoacrylate is purchased as a working solution. Shelf life is indefinite for this purchased product.

1. Place aluminum dish on a heating device and pour approximately one (1) teaspoon of glue in the dish. May use more or less glue depending on evidence being processed.
2. May add accelerator (ex: water), if required. Manufacturer recommends optimum humidity of 30-60%.
3. Place evidence into fuming chamber either by suspending or standing so all areas are exposed.
4. Seal fuming chamber.
5. Turn on heating device. Manufacturer recommends operation at 60°-85° F (16°-29° C) and not to heat above 250° F (120° C).
6. After latent print(s) are developed (usually 8 to 30 minutes), turn the heater off and exhaust the fumes from chamber before opening. May process item longer with superglue, as needed, for adequate clarity of results.
7. Vent chamber (usually 10 or more minutes).
8. Remove evidence and view for developed latent prints. Oblique and/or intense light may be utilized to better visualize developed latent prints.
9. If suitable latent prints are developed, the examiner may indicate the latent print with suitable markings as appropriate to be preserved by digital imaging/photography.
10. Depending on type of evidence, additional processing techniques for development of latent prints may be used (ex: Powder Processing and/or Fluorescent Dye Staining).

## Sequence

Cyanoacrylate Ester Fuming is typically utilized after a visual examination and before the utilization and application of other processes (ex: Powder Processing and/or Fluorescent Dye Staining). Cyanoacrylate Ester Fuming may interfere with DNA analysis and latent print blood enhancement techniques (such as Amido Black processing).

ARCHIVED 7/2/2012

## **DFO (1, 8 – DIAZAFLUOREN – 9 – ONE)**

### **Description of Process**

DFO is used to develop latent prints on porous surfaces. DFO is a Ninhydrin analogue and reacts to the amino acids present in perspiration. The prints will appear a pinkish-orange color; however, when viewed under various lasers and alternate light sources the prints will fluoresce brightly and are much more visible, especially on a dark colored surface that might hide prints that have been developed with Ninhydrin alone.

### **Safety**

Reagent is flammable. It is a sensitizer and causes staining of the skin. Mixing must be performed in a vent hood wearing lab coat, gloves, and eye protection. DFO is mixed with carriers that are highly flammable and irritant. Wear gloves, a lab coat, safety eyewear and use in a lab fume hood. Must be disposed of like any other flammable chemical.

### **Process for Use**

DFO crystals are purchased and the DFO Stock solution and DFO Working Solution are prepared in the laboratory.

1. Apply or dip the item in the DFO working solution for approximately ten (10) seconds, allow drying for approximately three (3) minutes.
2. Repeat the process.
3. Heat is then applied to the dried specimen by placing it in an oven that contains no humidity or use an iron with no steam. Heat for ten (10) minutes at 100 C (212 F).
4. View under a laser or other alternate light source as the developed prints may be invisible to the naked eye.
5. Examine item for latent prints and indicate the latent with suitable markings as appropriate to be preserved by digital imaging/photography.

### **Sequence**

If other processes are to be used on the same piece of evidence, DFO should be used after IND and before Ninhydrin or Physical Developer.



## IND (1,2-INDANEDIONE)

### Description of Process

1,2-Indanedione is a fluorescent amino acid reagent applied for developing latent prints on porous surfaces such as paper and cardboard. There is no discoloration or background staining evident on the 1,2-Indanedione processed samples that consistently appears when processing with DFO. 1,2-Indanedione may be used in place of DFO.

### Safety

Chemicals used in preparation and process are flammable and irritant. Avoid contact with skin and eyes. Wear proper protective equipment when preparing and processing items: lab coat, gloves, and safety glasses (goggles). Wear amber protective eye wear when viewing results under laser light.

### Process for Use

IND crystals are purchased and a working solution is prepared in the laboratory.

1. Apply the IND solution to an item by spraying, dipping or brushing.
2. Allow to dry for approximately three minutes.
3. After the IND has dried, place the processed item in a humidity chamber to accelerate the development process.
4. 10 minutes at 100 C and 60% relative humidity.
5. The best results obtained for the thermal paper samples were achieved by not accelerating the development and allowing them to develop naturally in the laboratory environment from 4 to 12 hours.
6. Developed prints are observed through an orange/amber viewing filter using a light source.
7. Examine the item for latent prints and indicate the latent with suitable markings as appropriate to be preserved by photography.
8. The results can also be seen on some samples with white light and develop as a light pale pink color.

### Sequence

If other processes are to be used on the same piece of evidence, IND should be used prior to DFO, Ninhydrin and Physical Developer.

# NINHYDRIN

## Description of Process

Ninhydrin is a chemical method for developing latent prints on porous surfaces and absorbent materials such as paper, cardboard, and smooth raw wood. This method is based on the reaction of Ninhydrin and amino acids that are present in latent print residue. The first known use of Ninhydrin for latent print processing was in the early 1950s. It is sensitive to old prints as well as fresh prints.

Ninhydrin can be mixed using two carriers: acetone or petroleum ether.

## Safety

Ninhydrin should be used in a laboratory fume hood, a well-ventilated area, or outside.

Gloves, lab coat, and safety eyewear must be worn when using.

Ninhydrin is mixed with a carrier such as methanol, acetone, or petroleum ether (which is highly flammable).

Excess is disposed of as any flammable liquid.

## Process for Use

Ninhydrin crystals are purchased, but a working solution is typically prepared in the laboratory. Shelf life is approximately six months for a working solution.

1. Select the appropriate Ninhydrin base solution dependent upon the other substances on the surface. Acetone will cause certain inks to dissolve. Therefore, handwriting analysis should be performed before processing for latent prints begins.
2. Apply Ninhydrin solution to an item by spraying, dipping, or brushing.
3. After the Ninhydrin has dried, place the processed item in humidity chamber or steam the item with an iron to accelerate the development process.
4. Examine the item for latent prints and indicate the latent with suitable markings as appropriate to be preserved by digital imaging/photography or scanning.

Note: Developed latents on some dark-colored surfaces may be viewed with a light source for increased contrast. Development of latent prints may vary with exposure time to Ninhydrin.

## Sequence

Evidence that may have potential DNA evidentiary value may be processed for latent prints with Ninhydrin previous to DNA sample collection.

Ninhydrin can be utilized by itself or in conjunction with other processes if used in the following order:

1. IND; 2. DFO; 3. Ninhydrin; 4. Physical Developer

# POWDERS - PLAIN

## Description of Process

Black powder is one of the oldest and most basic methods of developing latent prints on non-porous surfaces. Powder is applied to a surface by lightly dusting over the surface with a soft bristle type brush or duster. Once prints are developed they should be preserved for later comparison by either photography and/or lifted with lifting tape and placed on a lift card. Black powder is not the only color available; however, it is the most commonly used type of powder (even on dark colored surfaces). In addition, there are fluorescent powders that may be used, which require the use of an alternate light source and appropriate filters.

## Safety

Powder can easily be inhaled. Wear a facemask to filter out loose powder in the air or dust in a fume hood. Wearing gloves will prevent you from getting your hands dirty; however once lifting tape is placed over the developed latent; gloves make it difficult to rub out air bubbles.

## Process for Use

All powder is purchased and not prepared in the laboratory.

1. Prior to applying the powder to the surface that is to be processed, view the item for prints that are visible to the eye; such as prints in blood, grease, or any foreign residue. If there are visible latent, patent, or plastic prints:
  - a) *Do not handle the evidence excessively since the item has not yet been processed.*
  - b) *Consider alternative processing techniques for prints in blood, grease, or any foreign residue.*
2. Take the jar of powder and tap the jar into the palm of the hand several times to break up clogs of powder and loosen the powder that has settled.
3. Pour a small amount of powder (two to three tablespoons) into a container.
4. Choose a type of brush to apply powder.
  - a) *Fiberglass or nylon – to be used on small or large objects.*
  - b) *Feather duster – for larger objects.*
  - c) *Short bristle brush – for small objects and also used for cleaning up latents by lightly brushing in the direction of ridge flow.*
5. Hold fiberglass brush, nylon brush, or feather duster between palms of your hands, rub hands back and forth several times to loosen and fluff out bristles or feathers.
  - a) *Dip the brush into the container of powder lightly to pick the powder up.*
  - b) *Tap the brush several times with your index finger over container to release excess powder.*
6. Apply the powder to the surface by lightly dusting over the surface (only the tips of bristles or feathers should touch the surface).

*a) Twirling motion – fiberglass brush*

*b) Back and forth motion – fiberglass, feather duster, or short bristle brush*

*c) Figure eight type motion – feather duster*

*d) Once a latent is visible, view the latent and then apply a few more strokes of powder. If the latent starts to lighten up or starts looking spotty – stop processing. The latent is at its maximum contrast. Additional processing will destroy or deteriorate the latent.*

7. Remove the excess powder from the processed item.

*a) Tap item lightly on counter.*

*b) Use short bristle brush (brushing with the flow of ridges).*

*c) Make multiple lifts of the same print.*

8. Preserve the latent print for later comparison by either digital imaging/photography or a lifting technique appropriate to the evidence.

### **Sequence**

Powder Processing can be used at the Forensic Scientist or Technician's discretion. It should not be used on porous items such as checks, cardboard, paper sacks, etc. A chemical process would be best suited for these types of evidence; however, if the paper or cardboard has a shiny or slick surface (such as magazine covers or matchbook covers) it could be used.

## POWDERS – MAGNETIC

### Description of Process

With magnetic powder, there is no brush with fibers or bristles to hold the powder. The powder is actually made up of finely ground metal shavings with colored powder. The powder is applied with a metal rod or wand that has a magnet inside that attracts the powder-like whiskers. This method was developed in the early 1960's. It is not used for every piece of evidence; it is just another tool available in the latent processing field.

### Safety

There are no known safety hazards. Magnetic Powder Processing is not as messy as Powder Processing. Any overspill can be picked up with the wand and released back into the container.

### PROCESS FOR USE

All powder is purchased and not prepared in the laboratory.

1. Prior to applying the powder to the surface that is to be processed, view the item for prints that are visible to the eye; such as prints in blood, grease, or any foreign residue. If there are visible latent, patent, or plastic prints:
  - a) *Do not handle the evidence excessively since the item has not yet been processed.*
  - b) *Consider alternative processing techniques for prints in blood, grease, or any foreign residue.*
2. Take the jar of powder and tap the jar into the palm of the hand several times to break up clogs of powder and loosen powder that has settled.
3. Stick large bulb end of wand into the jar to pick up metal shavings. Pull the rod out of the wand to release metal shavings.
4. Go over the surface using a back and forth motion with only the metal shavings coming in to contact with the surface. CAUTION: If the metal bulb end comes in to contact with the surface, it could scratch or destroy a latent print.
5. Remove excess powder from the item processed.
  - a) *Tap the item lightly on the counter.*
  - b) *Use a short bristle brush (brushing with the flow of ridges).*
  - c) *Make multiple lifts of the same print as necessary.*
6. Preserve the latent print for later comparison by either photography or a lifting technique appropriate to the evidence.

### Sequence

Magnetic Powder can be used at the Forensic Scientist or Technician's discretion, usually after cyanoacrylate fuming. Magnetic Powder Processing is not suited for processing metal objects or porous items such as checks, raw cardboard, or paper sacks; however, it can be used on shiny slick surfaces such as magazine covers, match book covers, etc.

# PHYSICAL DEVELOPER PROCESSING

## Description of Process

The Physical Developer process is used for processing porous surfaces, especially on porous surfaces that have been wet, and on U.S. currency. It reacts with the fats, oils, and waxes present in the fingerprint residue. Until the introduction of a two-solution pre-mixed kit, Physical Developer had to be mixed from seven different chemicals following a complicated mixing routine.

## Safety

The Physical Developer process can be carried out with no known health hazards, provided a few precautions are carried out, such as wearing lab coats, latex gloves, and safety eyewear. The reagents in the working solutions are corrosive and toxic and will cause black staining on skin and clothing. Maleic acid is extremely destructive to tissue of the mucous membranes and upper respiratory tract, eyes, and skin. Mix in a fume hood while wearing gloves, lab coat and eye protection.

## Process for Use

Physical Developer is typically purchased as a working solution and not prepared in the laboratory.

1. Lay out three (3) glass trays
  - a) *Maleic acid prewash* If not used the Physical Developer will cause the paper to turn dark and obliterate any latent prints.
  - b) *Physical Developer working solution*
  - c) *Water rinse*
2. Immerse the item(s) in the prewash until the bubbles stop (do NOT use metal tongs).
3. Transfer the item(s) to the working solution. The tray with the working solution should be placed on an orbital shaker or manually rocked back and forth by hand. Leave the item(s) in the working solution for 5 to 15 minutes or until latents are developed.
4. Remove the item(s) and place them in a water rinse to remove excess solution.
5. Remove the item(s) from the water rinse and let dry.
6. Examine the item for latent prints and indicate the latent with suitable markings as appropriate to be preserved by photography.
7. If latents are hidden by a dark background color the background can be lightened up by placing the item(s) in a 50% household bleach – 50% tap water solution.
8. The prints will appear grayish-brown in color.

## Sequence

Physical Developer can be used in conjunction with other processes but normally after IND, DFO, and Ninhydrin.

# RHODAMINE 6G FLUORESCENT DYE PROCESSING

## Description of Process

Rhodamine 6G is to be used only on non-porous items and after the item has been treated with cyanoacrylate fuming. Rhodamine 6G is a fluorescent dye used to make cyanoacrylate developed latent prints more visible on various colored surfaces. Lasers or alternate light sources are used in conjunction with this process. Rhodamine 6G enhanced latent prints will have to be photographed under a light source. Different carriers for the working solution can be utilized to decrease processing times, preserve inked markings on evidence, or for use on special surfaces.

## Safety

Rhodamine 6G working solution and stock solutions are extremely flammable and caution should be used. This reagent should be mixed and applied to evidence under a fume hood so it is not inhaled. Gloves, a lab coat, and eye protection should be used.

## Process for Use

Rhodamine 6G crystals or powder is purchased and working and stock solutions are prepared in the laboratory.

1. After the evidence has been processed by cyanoacrylate ester fuming, apply the appropriate Rhodamine 6G working solution by either dipping, or using a spray device or squirt bottle.
2. Place the evidence under a fume hood to dry.
3. Examine the evidence under the laser and view using an orange filter. The power setting (beam intensity) may be adjusted as needed.
4. If the dye appears to be in excess, it may be rinsed with an application of distilled water or methanol over the evidence to reduce its thickness. A second application of dye stain may be necessary after the rinse.
5. Examine the item for latent prints and indicate the latent with suitable markings as appropriate to be preserved by photography.

## Sequence

Rhodamine 6G is to be applied only after a non-porous has been exposed to cyanoacrylate fuming. Powder application may be utilized before or after use of Rhodamine 6G.

# **WETWOP/WET POWDER**

## **Description of Process**

Wetwop and other sticky side powder equivalents are used to develop latent prints on the adhesive sides of tapes, decals, and other items.

## **Safety**

No known safety hazards. Although manufacturer caution states "Wet Powder is a highly stainable product. Use proper clothing."

## **Process for Use**

Wetwop and Wet Powder are typically purchased as a working solution. Shelf life is indefinite for purchased Wetwop and Wet Powder.

1. Shake container before use.
2. Pour small working amount of Wetwop, Wet Powder, or other sticky side powder equivalent into appropriately sized container.
3. Using latent print/fingerprint brush, apply solution onto the adhesive side of tape or other adhesive surface.
4. Leave on for 10-15 seconds.
5. Rinse off under slow running, cool water.
6. Let dry.
7. Examine the item for latent prints and indicate the latent with suitable detail as appropriate to be preserved by digital imaging/photography.

## **Sequence**

Wetwop/Wet Powder and other sticky side powder equivalents are typically used instead of other processes (cyanoacrylate ester fuming, Rhodmine 6G, etc.) on the adhesive side of items. The non-adhesive side may be processed as normal before applying Wetwop/Wet Powder and other sticky side powder equivalents. Adhesive surface processing may interfere with DNA analysis.



## **APPENDIX C - - WIN/AAFIS, IAFIS**

### **COMPUTER DATABASES**

WIN/AAFIS is the Western Identification Network/Alaska Automated Fingerprint Identification System. It is a ten print record database jointly shared by Alaska, Oregon, Idaho, Utah, Nevada, Montana, Washington, BICE (Bureau of Immigration and Customs Enforcement) and Wyoming. California and several other individual agencies are interface members of the network.

IAFIS is the Integrated Automated Fingerprint Identification System. It is the system and database maintained by the FBI.

The Latent Print Training Manual includes a section on the use of WIN. The procedures used can be found in the WIN AFIS21 Global Workstation GWS-L User Guide and Print Quest AFIS Manual.

IAFIS is not frequently utilized. Latent Examiners can be trained by a designated WIN Trainer. The procedures for IAFIS are outlined in the WIN AFIS21 Global Workstation GWS-L User Guide, Appendix B.

### **SEARCHING**

All latent prints that are of sufficient quality and have not been identified with known finger or palm prints can be entered into WIN/AAFIS, AFIS, or IAFIS if requested on the Request for Lab Services Form or if confirmed with the requesting agency/officer.

Latents are searched at the discretion of the examiner signing the report for the case.

### **INDIVIDUAL CHARACTERISTIC DATABASE SAMPLES**

The individual characteristic database sample files for the Latent Print Discipline contain known palm, major case prints, ten print cards and sole impressions from suspect and elimination persons. These impressions are reference samples.

The reference samples are stored in a locked room. Latent Discipline employees have access to the room. The files for each individual are stored in manila envelopes and arranged alphabetically by last name. The samples are uniquely identified by name and laboratory case number.

## APPENDIX D - - DNA SAMPLING

Analysts in the Latent Print Discipline may be responsible for collecting contact DNA from evidence items.

The analyst should rely on the DNA Section to determine if any of their processing substances may interfere with subsequent DNA analysis. The analyst may use their discretion to determine the order of processing, including at what stage the DNA sample is taken. This decision is based upon training and experience and is dependent upon the nature and the condition of the evidence.

Due to the possibility of DNA transfer which can occur with the re-use of fingerprint powder brushes, it is recommended that DNA be collected before the application of fingerprint powders.

### Isolating DNA Samples

1. Standard casework precautions should be observed to prevent sample contamination, such as clean work table, evidence on fresh sheet of paper, new gloves, mask and new and/or cleaned sampling materials.
2. A sample may be isolated by swabbing it.
3. Let the DNA analysts know about any chemicals or processing that was performed in this section on the sample area prior to sample collection. This information may be included on the DNA Case Status Summary form, and should appear in bench notes. It may also be noted on the sample packaging.

### Swabbing Method

1. Take two swabs of each sample area (eg. Mouth area of can/bottle, textured areas on firearm). The first swab applied to the stain should be wet, followed by a dry swab. These swabs should be packaged together as one sample.
2. Water used to moisten the swab should be prepared by the DNA section. Water used for swabbing should be replaced one month after being opened. The date opened is recorded on the bottle. This date and the lot number of the water should be in the analyst's case notes. The pipette system with the disposable pipette is the approved method for applying water to the swab for IN LAB swabbing.
3. If other water is used (not prepared by the DNA section), a control swab of the water should be prepared. This control swab should be packaged as a separate item.
4. Swabs should be air-dried before packaging.

## Packaging DNA Samples

1. Each separate area sampled should be given a unique number. For example, if there are two separate samples from item 41, then they should be labeled as 41-1 and 41-2 using hyphens to separate the item number from the sample number.
2. Each sample should be packaged separately and labeled.

Example for swab label:

**Lab # 00-00000**

**Swab 41-1 from the mouth of the bottle (Item 41)**

3. The packaged samples should then be placed in an outer envelope which is labeled with a complete description of the contents and labeled with the item number and the analysts initials (for example: as item #41LCH (LCH being the analyst's initials).
4. The created item should be added as evidence in LIMS and the DNA pending assignment in JTRAX is assigned to the DNA Discipline Supervisor to alert the DNA Discipline that contact DNA sampling has been completed.

## APPENDIX E - - REFERENCES

### Amido Black

Lee HC, and Gaensslen RE. Advances In Fingerprint Technology.

Federal Bureau of Identification, Latent Print Unit. Chemical Formulas and Processing Guide For Developing Latent Prints. 1999.

### Cyanoacrylate

Lee HC, and Gaensslen RE. Advances In Fingerprint Technology.

Wertheim PA. Atmospheric Super Glue Processing. Lightning Powder Minutiae. Sept-Oct, 1997.

### DFO

Lightning Powder. 1, 8 – Diazfluoren – 9 – One (DFO), No. 1-0038. Mar 1991.

### IND

Almog J, Springer E, Wiesner S, Frank A, et al. Latent Fingerprint Visualization by 1,2-Indanedione and Related Compounds: Preliminary Results. J For Sci. 1999. 44(1), pp. 114-118.

Hauze DB, Petrovskaia O, Taylor B, Joullie MM, Ramotowski R, and Cantu AA. 1,2-Indanediones: New Reagents for Visualizing the Amino Acid Components of Latent Prints. J For Sci. 1998. 43(4), pp. 744-747.

BVDA International. 1,2-IND, Catalogue No. B-78100. Online, [http://www.bvda.com/EN/prdctinf/en\\_ind\\_1.html](http://www.bvda.com/EN/prdctinf/en_ind_1.html)

Roux C, Jones N, Lennard C, and Stilovic M. Evaluation of 1,2-Indanedione and 5,6-dimethoxy-1,2-Indanedione for the Detection of Latent Fingerprints on Porous Surfaces. J For Sci. 2000. 45(4), pp 761-769.

Sirchie Finger Print Laboratories, Inc. 1,2-Indanedione, Catalog No. LV508. 2003.

Wiesner S, Springer E, Sasson Y, and Almog J. Chemical Development of Latent Fingerprints: 1,2-Indanedione Has Come of Age. J For Sci. 2001. 46(5), pp. 1082-1084.

### Ninhydrin

Lee HC, and Gaensslen RE. Advances In Fingerprint Technology.

Federal Bureau of Identification, Latent Print Unit. Chemical Formulas and Processing Guide For Developing Latent Prints. 1999.

### Physical Developer

Lightning Powder. Physical Developer Kit, No. 1-2730. Apr 2000.

## **Powders - Plain**

Olson RD. Scott's Fingerprints Mechanics.

Moenssens AA. Fingerprint Techniques.

Federal Bureau of Investigation. The Science of Fingerprints.

Wertheim PA. Black Powder Processing. Lightning Powder Minutiae. May-Jun, 1997.

## **Powders - Magnetic**

Moenssens AA. Fingerprint Techniques.

Wertheim PA. Magnetic Powder. Lightning Powder Minutiae. Jul-Aug, 1997.

## **Rhodamine**

Federal Bureau of Investigation, Latent Print Section. Processing Guide For Developing Latent Prints. 2000.

## **Wetwop**

Lightning Powder, Inc. Sticky-Side Powder, No. 1-2721.

Federal Bureau of Identification, Latent Print Section. Ash Gray Powder. Processing Guide For Developing Latent Prints. 2000.

Armor Forensics, Wetwop, No #1-0078.

Kjell Carlsson Innovations, Wet Powder.

## **WIN-AAFIS**

AFIS Latent Fingerprint Best Practices, Western Identification Network, Inc., October 2008

## APPENDIX F - - REVISION HISTORY

LP2011 R0 Effective date 08/24/11 - Total rewrite of manual for ISO standards

| R0<br>Page | Change made   | R0<br>Page | Change made  |
|------------|---|------------|--|
| 1          | Updated Document Structure to reflect Removal of <b>Appendix A</b> and new page numbers   | 4          | <u>Removed Abbreviations</u><br>EN Enhanced<br>FLP Fluorescent Powder<br>FPC Finger Print Card<br>GB Gun Blue<br>EN Enhanced<br>FLP Fluorescent Powder<br>ID2 Identified to<br>IO Iodine<br>IPP Inked Palm Print<br>LA Labeled As<br>LC Lift Card(s)<br>Lt Left<br><br>Replaced LtI with LI<br>Replaced LtL with LL<br>Replaced LtM with LM<br>Replaced LtR with LR<br>Replaced LtT with LT<br><br><u>Removed Abbreviations</u><br>LtP Left Palm<br>MAP Maleic Acid Prewash<br>MCP Major Case Prints |
| 3          | <u>Removed Abbreviations</u><br>Aba Amido Black- Aqueous<br>ABRa Amido Black Rinse- Aqueous<br>ABm Amido Black- Methanol Carrier<br>ABRm Amido Black Rinse- Methanol Carrier<br><b>Replaced with AB Amido Black</b><br><br><u>Removed Abbreviations</u><br>AQ AFIS Quality<br>BG Bag<br>BiP Bichromatic Powder<br>BiMP Bichromatic Magnetic Powder<br>BLK Black<br>BMP Black Magnetic Powder<br>Brn Brown<br>BY Basic Yellow<br>C: Containing<br>CHAR Characteristics<br>CHEM Chemistry<br>CLR Clear<br>COMP Compared or Comparison |            |  |
| 5          | <u>Removed Abbreviations</u><br>NINd Dry Ninhydrin Process<br>NINh Ninhydrin HFE Carrier<br>NINp Ninhydrin Petroleum Ether Carrier<br><b>Replaced with NIN Ninhydrin</b><br><br><u>Removed Abbreviations</u><br>P Palm<br>PG Plastic Bag<br>PH Photographed<br>Pi Pink<br>POSS Possible<br>PP Palm Print<br>PPC Palm Print Card<br>Rt Right<br><br>Replaced Rtl with RI<br>Replaced RtL with RL<br>Replaced RtM with RM<br>Replaced RtR with RR<br>Replaced RtT with RT<br><br><u>Removed Abbreviation</u><br>RiP Right Palm        | 6          | <u>Removed Abbreviations</u><br>SPRg Small Particle Reagent – Gray<br>SPRw Small Particle Reagent - White  |
|            |   | 10         | 4.6.2 Initial control testing of . . . chemicals mixed at the laboratory is noted in the JTRAX Chemical Inventory for each batch.<br><br>Changed to<br><br>4.6.2 Initial control testing of . . . chemicals mixed at the laboratory is noted in the CHEM INV Excel Spreadsheet for each batch.<br>-----<br>4.6.4 <b>Appendix A</b> lists the Latent Discipline Approved Vendors.<br><br>Changed to<br><br>4.6.4 The Latent Discipline does not have any critical consumables.                        |
|            |   | 12         | 4.13.1.2 Digital photos and/or scans.....flow chart is shown in <b>Appendix B</b><br><br>Changed to<br><br>4.13.1.2 Digital photos and/or scans.....flow chart is shown in <b>Appendix A</b>   |

| R0<br>Page | Change made | R0<br>Page | Change made |
|------------|-------------|------------|-------------|
|------------|-------------|------------|-------------|

|    |   |    |   |
|----|---|----|---|
| 19 | <p><u>4.13.2.12</u> All identifications must be verified</p> <p>Changed to</p> <p><u>4.13.2.12</u> All identifications and exclusions must be verified</p>  | 24 | <p><u>5.7.1</u> Usually, an analyst works all the evidence submitted for a case. However, the question being asked in the Latent Print Discipline is “Did an individual leave a latent print on an item of evidence?” An analyst can sample by quality or quantity and select the evidence that has the best chance of retaining latent prints or selecting and working a sampling of the evidence until the individuals of interest are identified. A technician, may sample evidence items if they consult with an analyst to compare the developed prints. Once an analyst identifies the individual(s) of interest on an item, processing and comparison can cease.</p> <p>Changed to</p> <p><u>5.7.1</u> The question being asked in the Latent Print Discipline is “Did an individual leave a latent print on an item of evidence?” An analyst can select samples by quality and/or quantity. An Analyst can select the evidence that has the best chance of retaining latent prints or an analyst can select and process a portion of the evidence until the individuals of interest are identified. A technician can employ sample selection if they consult with an analyst. Once an analyst identifies the individual(s) of interest on an item, processing and comparison can cease.</p> <p>-----</p> <p><u>5.7.2</u> Unusual sampling plan situations must be approved by the Latent Print Supervisor.</p> <p>Changed to</p> <p><u>5.7.2</u> Unusual sample selection situations must be approved by the latent discipline supervisor.</p> <p>-----</p> <p><u>5.7.3</u> An analyst or laboratory technician.....at a minimum.</p> <p>Changed to</p> <p><u>5.7.3</u> An analyst or laboratory technician must document sampled selection in their notes. When sample selection is employed, the selected items must be identifiable at a later date. An analyst or laboratory technician must initial the sample selection items in a case at a minimum. If this is not possible, the sampled selection items may be repackaged and the analyst or technician must, at a minimum, initial the packaging.</p> <p>-----</p> <p><u>5.8.4.2</u> Evidence not in the process of examination is stored a locked tote. Larger evidence may be stored in the latent case archive room.</p> <p>Changed to</p> <p><u>5.8.4.2</u> Evidence not in the process of examination is stored in a locked tote or evidence cabinet. Larger evidence may be stored in the latent case archive room.</p> |
| 20 | <p><u>5.1.3.1</u> Initial control testing of Rhodamine . . . . chemicals mixed at the laboratory is noted in the JTRAX Chemical Inventory for each batch.</p> <p>Changed to</p> <p><u>5.1.3.1</u> Initial control testing of Rhodamine. . . . chemicals mixed at the laboratory is noted in the CHEM INV Excel Spreadsheet for each batch.</p>  |    |   |
| 21 | <p><u>5.4.1</u> Mixing instructions for reagents made at the laboratory are in the JTRAX Chemical Inventory.</p> <p>Changed to</p> <p><u>5.4.1</u> Mixing instructions for reagents made at the laboratory are in the CHEM INV Excel Spreadsheet.</p> <p>-----</p> <p><u>5.4.2</u> Processing used for a case evidence ..... are listed in <b>Appendix C</b> of this manual.</p> <p>Changed to</p> <p><u>5.4.2</u> Processing used for a case evidence.....are listed in <b>Appendix B</b> of this manual</p> |    |   |
| 22 | <p><u>5.5.3</u> Equipment manuals are stored in the Equipment folder in the Latents Share folder.</p> <p>Changed to</p> <p><u>5.5.3</u> Equipment manuals are stored in the Latents Share folder.</p> <p>-----</p> <p><u>5.5.5</u> A spread sheet with equipment records is stored in the Equipment folder in the Latents Share folder.</p> <p>Changed to</p> <p><u>5.5.5</u> Spread sheets with equipment records is stored in the Latents Share folder.</p>   |    |   |

| R0<br>Page | Change made   |
|------------|---|
| 25         | <u>5.8.4.6.1</u> <i>Nothing Additional</i><br>Changed to<br><u>5.8.4.6.1</u> Palm files are treated as Reference Material                     |
| 27         | Removed<br><b>APPENDIX A - - APPROVED VENDORS</b>   |
| 28         | Renamed<br><b>APPENDIX B - - PHYSICAL SECTION IMAGES FOLDER FLOW CHART</b><br><b>APPENDIX A - - PHYSICAL SECTION IMAGES FOLDER FLOW CHART</b> |
| 29         | Renamed<br><b>APPENDIX C - - PROCESSING WORK INSTRUCTIONS</b><br><b>APPENDIX B - - PROCESSING WORK INSTRUCTIONS</b>                           |
| 43         | Renamed<br><b>APPENDIX D - - WIN/AAFIS, IAFIS</b><br><b>APPENDIX C - - WIN/AAFIS, IAFIS</b>   |
| 44         | Renamed<br><b>APPENDIX E - - DNA SAMPLING</b><br><b>APPENDIX D - - DNA SAMPLING</b>   |
| 46         | Renamed<br><b>APPENDIX F - - REFERENCES</b><br><b>APPENDIX E - - REFERENCES</b>   |
| 43         | Renamed<br><b>APPENDIX D - - WIN/AAFIS, IAFIS</b><br><b>APPENDIX C - - WIN/AAFIS, IAFIS</b>   |
| 48         | Renamed<br><b>APPENDIX G - - REVISION HISTORY</b><br><b>APPENDIX F - - REVISION HISTORY</b>   |